

Autoclavable Micro Pipette

Nichipet Premium

5-year warranty


Innovation in durability, strength, stability and comfort the latest model in Nichipet series



Enlarged and easy to read counter display for volume setting.




Non corrosive ceramic used in the plunger. (100 μ L - 10000 μ L)

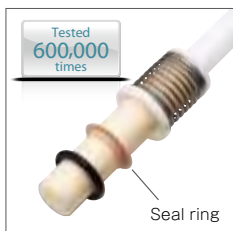


Seal ring for the Premium is tested dispensing / aspirating 600,000 times. Non-greased plunger makes it easy to maintain. (2 μ L - 1000 μ L)
The test was conducted by Nichiryo. The test result may vary with condition of usage.



New material used for 4 notch-lock improves its strength, which minimizes the lever shift from the original position.

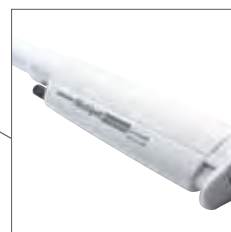

Materials used for pipette body have high tolerance to solvents and drop impact.

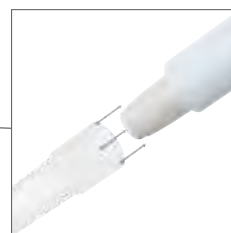

Nozzle top has highly abrasion tolerant. (2 μ L - 1000 μ L)

Innovation in "Durability" and "Strength"

- Seal ring for the Premium is tested dispensing / aspirating 600,000 times. (2 μ L - 1000 μ L)
- Nozzle top is highly abrasion tolerant. (2 μ L - 1000 μ L)
- Non corrosive ceramic is used for the plunger. (100 μ L - 10000 μ L)

Innovation in "Stability" and "Preciseness"

- Hyper Blower System, longer second-push stroke, improves dispensing efficiency. (2 μ L, 10 μ L, 20 μ L)
- Easy-Calibration function provides simple lab calibration.
- Hand temperature doesn't permeate through the body of the pipette which works to avoid affecting volume measurement. (Patented)

Innovation in "Fatigue-Free" and "Comfort"


- New round shape of the Premium mitigates fatigue and stress from heavy workload.
- Enlarged push button makes it easy to depress.
- White body color mitigates reflection of room light.

*All features of Nichipet EXII are available on Nichipet Premium.

What's new!

Shortened to release tension to your thumb.

Shortened the distance to operate smoothly.

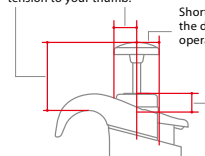

Specifications

Cat. No	Volume Setting	Volume range(μ L)	Increments(μ L)	Usable tips
00-NPP-2	Digital Setting	0.1 – 2	0.002	BMT2-UT/UTWR/UTWRB, FUT/FUTRB
00-NPP-10		0.5 – 10	0.01	BMT2-SS/SSR/SSRB, FSS/FSSRB
00-NPP-20		2 – 20	0.02	BMT2-SG/SGR/SGRB, FSG/FSGRB
00-NPP-100		10 – 100	0.1	BMT2-SG/SGR/SGRB, FSE/FSERB
00-NPP-200		20 – 200	0.2	BMT2-SG/SGR/SGRB, FAG/FAGRB
00-NPP-1000		100 – 1000	1	BMT2-LG/LGR/LGRB, FLG/FLGRB
00-NPP-5000		1000 – 5000	10	BMT2-X/XR/XRB
00-NPP-10000		1000 – 10000	10	BMT2-Z/ZR/ZRB

* See p.15-16 for tip selection.

Accuracy and Precision

Cat. No	Volume(μ L)	Accuracy(%)	Precision(%)	Cat. No	Volume(μ L)	Accuracy(%)	Precision(%)
00-NPP-2	0.2	$\pm 12.0^*$	$\leq 6.0^*$	00-NPP-200	20	± 1.0	≤ 0.5
	1	± 5.0	≤ 2.5		100	± 0.8	≤ 0.3
	2	± 3.0	≤ 1.0		200	± 0.8	≤ 0.2
00-NPP-10	1	± 4.0	≤ 3.0	00-NPP-1000	100	± 1.0	≤ 0.5
	5	± 1.0	≤ 1.0		500	± 0.8	≤ 0.3
	10	± 1.0	≤ 0.5		1000	± 0.7	≤ 0.2
00-NPP-20	2	± 5.0	≤ 3.0	00-NPP-5000	1000	± 1.0	≤ 0.3
	10	± 1.0	≤ 1.0		2500	± 0.8	≤ 0.3
	20	± 1.0	≤ 0.4		5000	± 0.6	≤ 0.2
00-NPP-100	10	± 2.0	≤ 1.0	00-NPP-10000	1000	± 2.0	≤ 0.4
	50	± 1.0	≤ 0.3		5000	± 0.8	≤ 0.3
	100	± 0.8	≤ 0.3		10000	± 0.4	≤ 0.2

*The AC and CV above are measured with the disposable tips described in the catalogue.

*The AC and CV at 0.2 μ L depend much on the sampling skill and the sampling condition where pipettes are used.

Nichiryo's quality starter packages

Nichipet Premium S-PAC

- Package type of 3 selected Pipettes with racked-Tips.
- Ready-to-start packages with 3 sizes of pipettes to perform high AC and CV pipetting.
- 3 packages are available to cover wide range of volume.



Nichipet Premium

3 pcs

Racked Tip

3 or 4 racks

Tube Opener

3 pcs



Nichiryo Tube Opener (code:00-NTO)

“Easier” and “safer” to open tube lid with holding pipette



- Attached on grip end, 0.5, 1.5 mL tube and screw cap can be opened easily.

*Compatible with Nichipet premium (2-1000 μ L) and Nichipet EXII (2-1000 μ L)

Nichipet Premium S-PAC

Description	Cat. No	Content
Nichipet Premium S-PAC A	00-SP-AP	NPP-10, NPP-100, NPP-1000 BMT2-SSR(0.5~10 μ L), BMT2-SGR(2~200 μ L), BMT2-LGR(100~1000 μ L), NTO
Nichipet Premium S-PAC B	00-SP-BP	NPP-20, NPP-200, NPP-1000 BMT2-SGR(2~200 μ L) x2, BMT2-LGR(100~1000 μ L), NTO
Nichipet Premium S-PAC H	00-SP-HP	NPP-2, NPP-10, NPP-100 BMT2-UTR(0.1~10 μ L) x2, BMT2-SSR(0.5~10 μ L), BMT2-SGR(2~200 μ L), NTO